CASADATM
healthcare

< relaxmyway

SCOOP

< relaxmyway

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Congratulations!

With the purchase of this massage device you have shown awareness of your health.

In order for you to enjoy the advantages of your massage device in the long term, we ask you to read and adhere to the instruction manual carefully.

We wish you a lot of fun with your personal Scoop.

Table of contents

Safety instructions	14
Equipment / Functions.....	16
Features.....	17
Technical Data	17
Scope of delivery	17
Using method	18
Maintenance instructions.....	20
Storage	20
Troubleshooting.....	20
Conformity declaration	21
Warranty conditions	21

< relaxmyway

Since being established in 2000, the name CASADA stands for high quality products on the health-spa and fitness market. Today, CASADA-equipment is available in 37 countries world-wide.

Constant further development, an unmistakable design and a focused orientation towards a better quality of life are what characterize CASADA-products, as much as their unbeatable and reliable quality.

Safety instructions

Please carefully read the following instructions before using this massage device to ensure its proper function and optimal effect. Please retain these instructions for use!

- Please use this product according to this manual
- Please ensure that the mains power is suitable for this product, do not use this product if its casing or mains cable is damaged
- Do not allow any foreign matter to enter the clearance between the massage head and the casing
- Do not allow the product to come in contact with water so as to avoid any possible accidents that may result from a short circuit
- Please disconnect the product from the mains supply after use or before cleaning to avoid damage to the equipment or personal injury
- Do not plug or unplug this product when your hands are wet
- Do not pull the cable when plugging or unplugging and don't handle it violently
- Do not damage the cable or change the circuit of the product
- Do not use any spare parts and attachments other than those recommended
- Do not use this product with too much force so as to avoid injury
- Do not use this product within an hour after a meal
- Please stop using this product and consult your local supplier if you notice anything abnormal with the product during use
- This appliance is not intended for use by persons (including children under the age of 14) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- If the power adapter is damaged, it must be replaced by the manufacturer, the maintain department or similar professional person in order to avoid danger.


Safety instructions

I. The following persons should only use the device after consulting a doctor:

- Persons who are mildly affected by electro-magnetic waves, such as users of pace makers or similar devices.
- Persons suffering from a malignant tumour.
- Persons with heart disease.
- Women who are menstruating, are pregnant or have just given birth.
- Persons suffering from diabetes and thus reduced sensory functions.
- Persons with skin lesions in the part of the body to be massaged.
- Persons with a body temperature of more than 38°C (fever).
- Persons suffering from acute inflammation / injuries at advanced stages (with symptoms such as tiredness, shivering, labile blood pressure, etc.).
- Persons suffering from extreme exhaustion.
- Persons suffering from osteoporosis, bruises or pulled muscles, fractures etc.
- Persons who are sensitive to heat and pain.
- Persons who cannot consciously move their legs and feet.
- Persons who are not capable of expressing their thoughts clearly and persons who have no or limited movements should not use this device.
- Persons who have been forbidden from having massages by the doctor should refrain from using this device.
- Persons with thrombosis, severe aneurysm, acute venous aneurysm, all types of skin rashes and skin diseases, etc.

Equipment / Functions

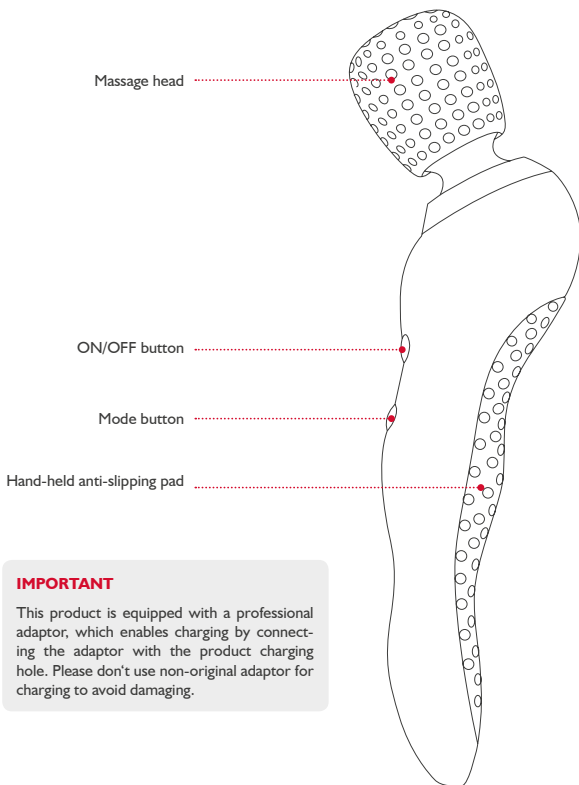IMPORTANT

This product is equipped with a professional adaptor, which enables charging by connecting the adaptor with the product charging hole. Please don't use non-original adaptor for charging to avoid damaging.

Features

Specials

 Vibration

 Automatic

 Wireless

Massage areas

Neck, back, waist and other suitable massage parts of human body

Special functions

Portable without wire
Travel lock
Flexible massage head

Technical data


Massage device	SCOOP
Product measurements	272 x 59 x 89 mm
Product weight	0,35 kg
Voltage adapter	
Input	AC 100-240 V ~ 50/60 Hz
Output	DC 5.0 V, 1500 mA
Voltage device	Lithium battery 3.7 V, 2200 mAh
Rated power	6 W
Automatic runtime	5 minutes
Certificates	   

Scope of delivery

- Massage device
- Cable with USB plug

Using method

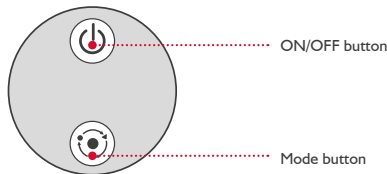ON/OFF BUTTON


When short-pressing this button the massage device turns on. When the device is turned on, light turns blue and the massage device runs at low speed. To turn off the massage device do the same combination. The blue light turns off and your device switches off.

MODE BUTTON

When the device is turned on, light flashes and turns blue by short pressing. There are 5 modes that can be chosen as shown on the next page. The blue light will flash shortly to confirm the mode change.

CHARGING

Insert the adaptor plug into the power socket, the other head of the plug into the product's charging hole and the light turns red; when the device is finished charging, the light turns green and you can use the product by pressing the ON/OFF button after unplugging the adaptor.

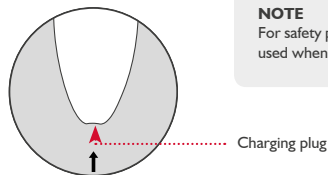NOTE

For safety purpose, the product cannot be used when charging

Using method

KEYLOCK

To avoid unintended operation or losing power by pressing the buttons, this product is designed with a travel lock; any operation under this mode is not possible. To unlock, press ON/OFF and mode button at the same time for longer than 3 seconds until blue light flashes once and vibrates once. This key combination acts both ways – locking and unlocking

NOTE

The product is key locked after purchase. Please unlock the massager before use.

Automatic programme	Speed / Vibration
Mode 1	<p>A graph with a vertical axis labeled '20%' and '0%'. A horizontal line is drawn at the 20% level.</p>
Mode 2	<p>A graph with a vertical axis labeled '50%' and '0%'. A horizontal line is drawn at the 50% level.</p>
Mode 3	<p>A graph with a vertical axis labeled '100%' and '0%'. A horizontal line is drawn at the 100% level.</p>
Mode 4	<p>A graph showing a series of pulses. The pulses start with low amplitude, gradually increase to a peak, and then gradually decrease in amplitude.</p>
Mode 5	<p>A graph showing a series of pulses with varying amplitudes and frequencies, creating a complex, irregular waveform.</p>

Maintenance instructions

Clean this product with neutral detergent and soft cotton. Gently rub the surface and wipe with a dry cloth. Please make sure the cloth is not too wet and don't use sharp brushes and detergents.

Storage

Store this massage device in a dry and cool area. Avoid excessive heating by direct sunlight as well as a place with high humidity.

Troubleshooting

Defect phenomenon	Possible causes	Solutions
Can't be turned on.	Low power or product is under travel lock.	a) Use after charging. b) Unlock the travel lock according to operation instruction.
Product shuts off automatically while using.	a) Battery power runs low. b) Auto timing (5 minutes) is over.	a) Use again after charging. b) Restart.
Product makes noise while using.	Normal sounds caused by motor.	This noise is normal and does not indicate a malfunction.

Conformity declaration

Corresponds to the European guidelines EMC 2014/30/EU, LVD 2014/35/EU.

Warranty conditions

The dealer grants the statutory warranty period on the product described overleaf. The warranty applies from the date of purchase. The date of purchase is demonstrated by the sales receipt.

The dealer is to repair or replace products registered within Germany which are found to be defective free of charge. Parts subject to wear and tear such as cloth covers are not covered. To assert your claim, simply return the defective product with the sales receipt before the end of the warranty period.

The warranty claim does not apply if it is established that a defect has occurred e.g. through external influence or as a result of repair or alteration not carried out by the manufacturer or an authorized dealer.

The warranty granted by the seller is limited to repair or replacement of the product. Under this warranty, the manufacturer or seller has no further liability and is not responsible for damage caused by disregard of the operating instructions and/or improper use of the product.

If inspection of the product by the seller reveals that the warranty claim relates to a defect not covered by the warranty or that the warranty period has expired, the cost of the inspection and repair are to be borne by the customer.

Casada International GmbH
Obermeiers Feld 3
33104 Paderborn, Germany
info@casada.com | www.casada.com

Copyright Images and texts are subject to the copyright of Casada International GmbH and may not be used for other purposes without express confirmation.

Copyright ©2017 Casada International GmbH. All rights reserved.

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Copyright ©2017 CASADA International GmbH. Alle Rechte vorbehalten. Alle Rechte vorbehalten.
Copyright ©2017 CASADA International GmbH. Tous droits réservés.
Copyright © 2017 Casada International GmbH. Wszystkie prawa zastrzeżone.